dyras

BTPS-2401B

HU	3
GB	8
DE.....	13
CZ.....	18
HR	23
SK.....	28
SI	33

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***d*yras**

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Modell: BTPS-2401B

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Personal scale
Model: BTPS-2401B

With „OKOK international” app

Please carefully read the instructions before using!

IMPORTANT SAFETY INFORMATION

You must obey the safety instructions that follow when using the appliance. Please read this section very carefully!

Of course, it is not possible to create safety regulations for every life situation, which is why we ask you to take all reasonable precautions in any given situation, even when you are taking the utmost care, in order to ensure safe operating conditions.


- This device is not supposed to be used by physically or mentally handicapped people, or persons without experience or knowledge (including children), unless their use of the device is supervised or guided by someone who is at the same time responsible for their safety. As for children, supervision is recommended to make sure that they do not play with the device. This device is not supposed to be used by children under the age of 14, not even with supervision!
- Never use an appliance or accessory that the manufacturer did not supply with, or propose for use with, the appliance. The use of such appliances or accessories might lead to accidents or irregular function.
- Do not place the appliance near hot objects (e.g. electric or gas cooker, hob, radiator, etc.) or on a hot surface and do not use the appliance on an unstable surface, or a high surface, or any surface where it is unstable and can be knocked down in an accident.
- Do not use the appliance in the presence of explosive or inflammable gases, vapours or liquids.

CAUTION

- Handle the balance with care, as it is a precise measuring instrument. Do not drop it or jump on it.
- Do not load the appliance with more than 180 kg to avoid damaging the device.
- Do not leave your scale with water on it, it may damage the internal parts.
- If the scales are not going to be used, please store them only in a horizontal position.
- Use your scale with caution. This is a high precision strain gauge sensors system.
- Do not get on it with wet feet, you could fall down
- The body fat analysis scale is not suitable for the following persons:
 - Children under 10 years and adults over 99 years
 - Persons with symptoms of fever, edema, or osteoporosis
 - Persons undergoing treatment for dialysis
 - Persons who use a pacemaker
 - Persons who take cardiovascular medication
 - Pregnant women
 - Sportspersons who perform over 5 hours of intensive sporting activities each day
 - Persons with a heartbeat lower than 60 beats/minute

OKOK INTERNATIONAL APP

Support Android 4.3 and above for selected Android devices with Bluetooth 4.0 and above. Install the app. Connect the unit with your mobile or any other device via Bluetooth. The analysis will be shown in the app called "OKOK" which you can easily download from the app store or google play. The app can save data for up to 8 people, so you and your family and friends can use this scale.

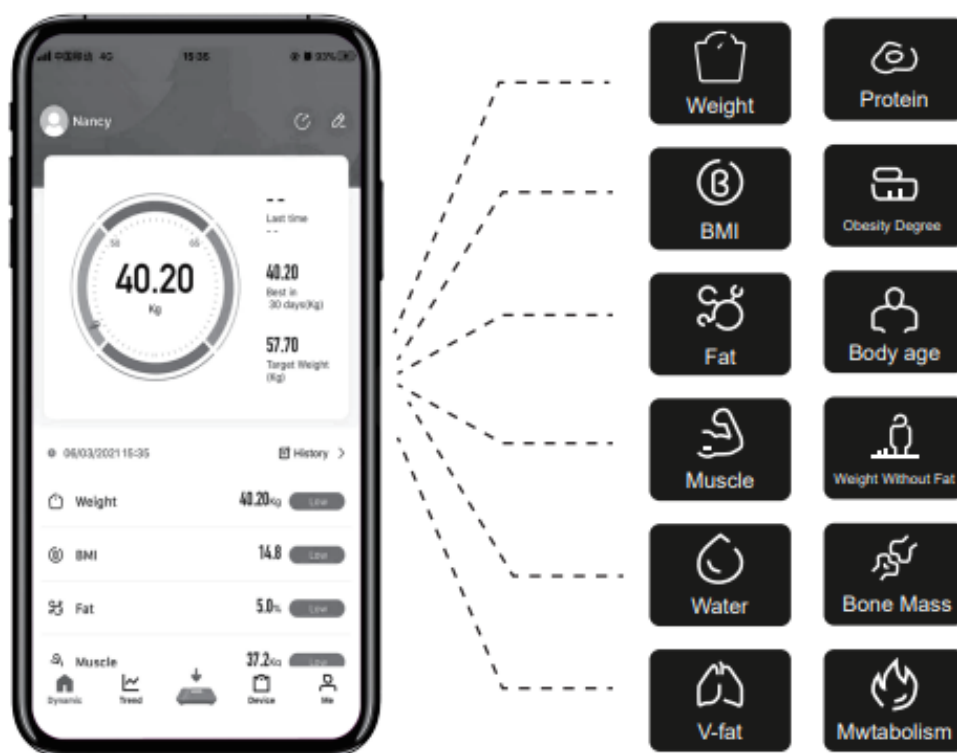HOW TO USE THE SCALE

The device needs to be connected with the app before first time use. The weight scale will automatically connect to your phone and the application. Make sure the Bluetooth is activated on your phone before you start using the weight scale.

1. Switch on Bluetooth Use your device to scan QR Code or search OKOK International from App Store/ Google Play to install APP.
2. Pair the device, register and complete user information (insert name, gender, height, birthday and preferred weight). Click next to confirm your personal information.

3. Weigh and synchronize. Step on the scale and weight yourself. Please wait till the device calculates your measurements.
4. View the result. The device will display body weight, BMI, fat, muscle, water, visceral fat, bones mass and metabolism. The device will save your details automatically.

Open the app and connect to Bluetooth. Pair the device.

- The user can type weight data into the app to keep a track.
- The device will display current weight and generate reports and give health advise based on measurements. Including current weight interface, weight list interface, date report and notice.
- Weight trend will include weekly report, monthly report and quarterly report, indicating changes of weight through curves.
- Account management includes registration, login, changing password, forget password, links to social media, and more.
- User management includes user information like nickname, birthday, height, gender, etc.

INSTRUCTIONS

- Remove the insulation film from the back of the scale near the battery compartment prior to operation.
- Place the batteries in the scale.
- Place the scale on a plane and hard surface prior to weighting.
- Pair your device with scale.
- Step on the scale, balance your weight evenly and wait until the display reads your weight. Your weight will flash on the readout for 10 seconds.
- Press the button near the battery box and you can choose the different unit: kg/lb

SPECIFICATIONS

- Power: 2XAAA battery
- Capacity: 180kg/396lb
- Unit: kg/lb
- Memory for 8 persons
- Measure weight: BMI/fat/V-fat/water/muscle/bone/Calorie
- Turn on auto
- Auto off in 10 seconds
- Low battery display: "Lo"
- Overweight display: "oL"

WASTE MANAGEMENT

Correct Disposal of this product	
	<p>This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.</p>

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CE SYMBOL

The mark is intended to indicate that the product complies with the relevant European Union standards and can be freely marketed in the European Union.

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Copyright © 2023 The dyras wwe.factory LLC U.S.A.
®= registered trademark of Dyras Europe Kft.