## RelaxGun G1

---

**EN – USER MANUAL**

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**HR – UPUTE ZA UPOTREBU**

---


### PACKAGE CONTENTS:

- 1× TrueLife RelaxGun G1
- 4× Massage head
- 1× USB-C cable
- 1× Manual
- 1× Safety instructions


### DESCRIPTION OF MESSAGE GUN:

1. Interchangeable massage head
2. Air outlet
3. USB charging port
4. Control panel
5. Handle

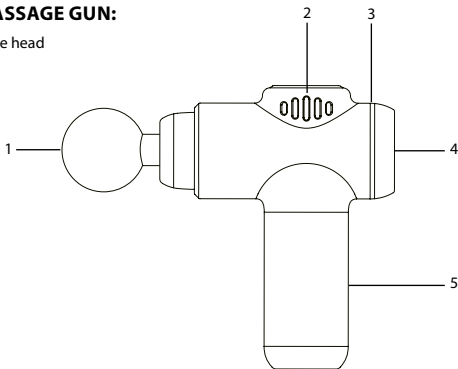
### DESCRIPTION OF CONTROL PANEL:

1. Display
2. – button (reduce intensity)
3. + button (increase intensity)
4. Switch

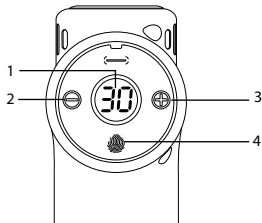
## DESCRIPTION AND USE OF MASSAGE HEADS:

### ROUND HEAD

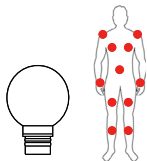
- Suitable for massaging almost all muscle areas on the body.
- The round shape of the head easily adapts to the shape of the body and can be used for massages targeted at large surface muscle groups.
- Suitable for use before and after exercise and at any time during the day.
- The head helps to gently relax muscle tissues, regenerating and stimulating individual points.
- Suitable intensity 1–30.

### FLAT HEAD

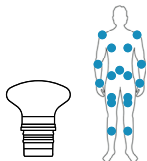
- Suitable for massaging almost all muscle areas on the body.
- The flat design and even pressure make it ideal for stimulating muscles and soft tissues near bones and joints.
- Suitable for use before and after exercise and at any time during the day.
- The head helps to gently relax muscle tissues, and regenerate and stimulate individual points including soft tissues in the joint area.
- Suitable intensity 1–30.

### FORKED HEAD

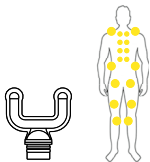- Suitable for indirect massage and stimulation of secondary points, adjacent muscles and deep muscle layers such as the trapezius muscles, Achilles tendons, calves, buttocks, cervical spine, ribs and cervical vertebrae.
- Suitable for use before and after exercise and at any time during the day.
- The head helps to gently relax muscle tissues, regenerating and stimulating secondary and deep points.
- Suitable intensity 1–15.

### POINTED HEAD

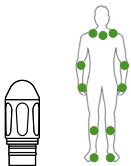- Suitable for direct point massage, deep stimulation of muscle layers and for parts of the body such as the arch of the feet, wrists, palms, back muscles, trapezius, neck muscles, forearms and calves.
- The head is suitable for use before and after exercise and at any time during the day.
- The head helps with targeted massage of muscles and soft tissues, regenerating and stimulating hard-to-reach points.
- Suitable intensity 1–10.

## **CHANGING THE HEAD**

Grasp the body of the gun with one hand and the head with the other hand.

Pull the head away from the body of the gun.

To put the head on, insert the end of the head into the hole in the gun. Gently push until you feel resistance.

## **SETTING THE INTENSITY:**

1. Hold the switch for two seconds to turn on.
2. Use the +/- buttons to adjust the intensity from 1 to 30, with 1 being the lowest intensity and 30 the highest.
3. When finished massaging, hold the switch for two seconds to turn off.

## **MESSAGE:**

1. Select and fit the appropriate massage head following the instructions in the "**DESCRIPTION AND USE OF MESSAGE HEADS**" section.
2. Hold the switch to turn on the massage gun.
3. Set the desired intensity. We recommend starting the massage with a lower intensity and gradually increasing it.
4. Apply gentle pressure to start massaging the desired part of the body.
5. Perform the massage by moving with the blood circulation, i.e. in the direction from the heart to the limbs.
6. The length of the massage should not exceed 15 minutes to avoid excessive strain on the muscles, bones and joints.

### **Note:**

If you experience any health complications, consult your doctor about the use of the product.

Do not use the product on the parts where you feel pain or have a skin injury.

Do not use the product on the face, spine or genitals.

Choose the appropriate massage head to suit your needs and health to maximise the benefits and minimise the risks.

Excessive pressure may cause injury to muscles or joints.

### **Warning:**

Do not use the product while it is charging.

The product will automatically stop after running continuously for 20 minutes. Press the + button to restart it.

### **CHARGING:**

Only charge the massage gun using the supplied USB-C cable.

Charge the massage gun before using it for the first time. Charging usually takes 4 hours.

Insert the end of the charging cable into the charging port on the massage gun. Plug the other end of the charging cable into a suitable mains adapter (recommended voltage 5 V / 2 A, adapter not included).

## SPECIFICATIONS:

Amplitude	6 mm
Number of intensity levels	1–30
Movement frequency	1,800–3,200 rpm
Maximum pressure	6 kg
Batteries	1500 mAh (7.4 V, lithium)
Battery life	3 hours
Full battery charge	3–4 hours
Input voltage	5 V, 2 A
Noise level	< 55 dB

## SAFETY INSTRUCTIONS:

Please read the enclosed safety instructions carefully.

For proper handling of the product, download and read the most up-to-date version of the manual, which can be downloaded at: [www.truelife.eu](http://www.truelife.eu).

**Printing errors and changes to information reserved.**

## EU DECLARATION OF CONFORMITY:

elem6 s.r.o. hereby declares that the Truelife RelaxGun G1 device is in compliance with the basic requirements and other relevant established directives RoHS 2011/65/EU Annex II 2015/863 and 2014/30/EU.

All Truelife products are intended for sale without restriction in Germany, the Czech Republic, Slovakia, Poland, Hungary and other EU member states. The full version of the Declaration of Conformity can be downloaded at [www.truelife.eu/support/doc](http://www.truelife.eu/support/doc).

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